


Gorilla[®] Breakaway Screw System

PRODUCT INFORMATION

Designed to Provide a Simple Solution for Syndesmosis Repair

- Once a patient begins weightbearing, traditional syndesmotic screws may fracture.
- The location of screw fracture of a traditional syndesmotic screw can vary, leading to difficult removal situations and/or potential pain at the screw fracture site.
- The Gorilla[®] Breakaway Screw System was designed so that if screw fracture occurs, the screw breaks cleanly at the notch point in the clear space.

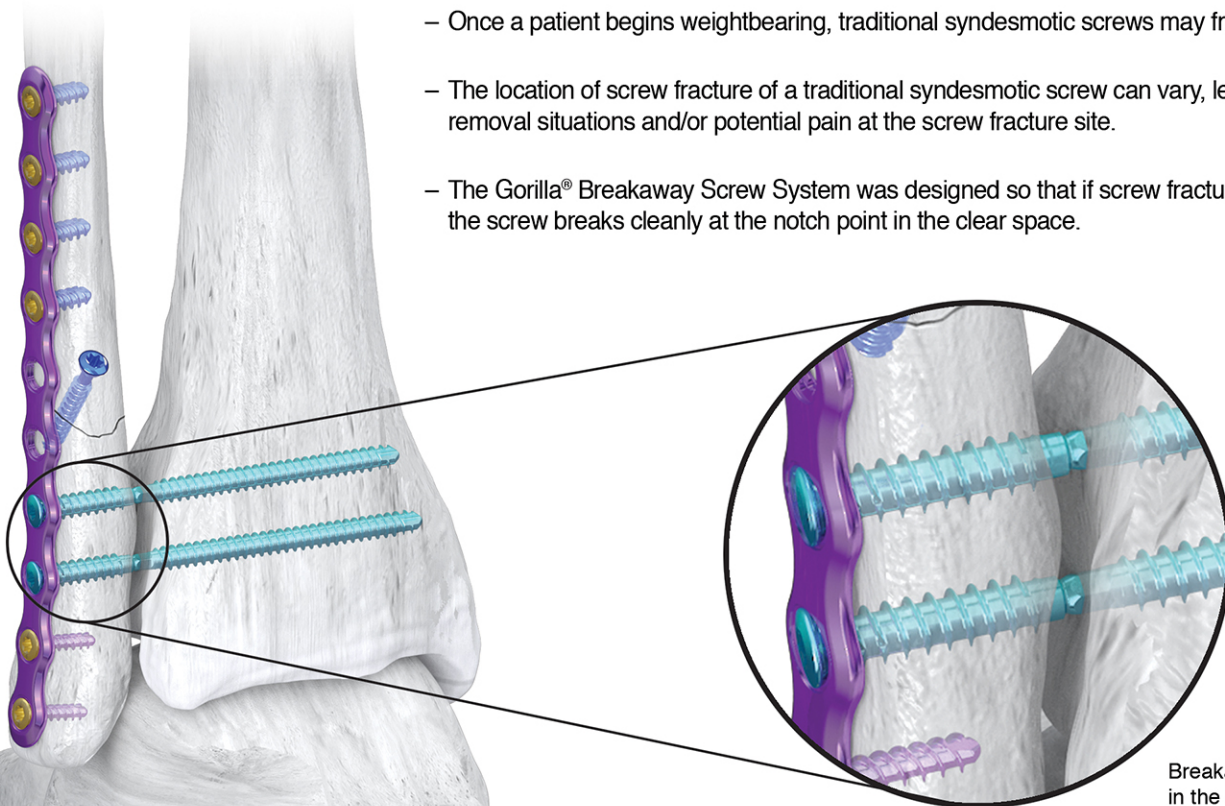

Breakaway notch resides
in the clear space between
fibula and tibia


Gorilla® Breakaway Screw System

IMPLANT OFFERING

BREAKAWAY™ SCREW SYSTEM

BREAKAWAY SCREW SIZES


3.9 mm Breakaway Screw - 14 mm Notch Length

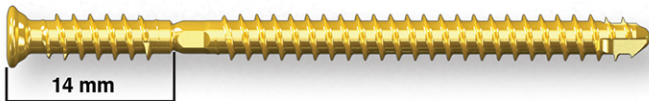

Screw Lengths:

40 mm - 50 mm (in 5 mm increments)
52 mm - 66 mm (in 2 mm increments)

3.9 mm Breakaway Screw - 17 mm Notch Length

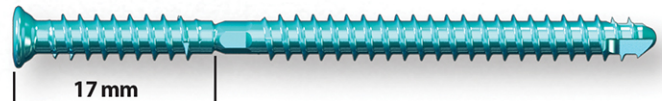


Screw Lengths:

40 mm - 50 mm (in 5 mm increments)
52 mm - 66 mm (in 2 mm increments)

SCREW FEATURES

- **Multiple Notch Lengths: 14 mm and 17 mm**
 - Addresses differences in patient anatomy and accommodates use with or without a plate
 - 2 mm increments from 52 - 66 mm intended for quadricortical fixation
- **Strength in Design**
 - Static bending strength and stiffness of the Gorilla® Breakaway Screw exceeded that of a comparable 3.5 mm solid screw as well as a 4.0 mm cannulated screw made of the same material.¹
- **If removal is desired, multiple removal features have been incorporated:**

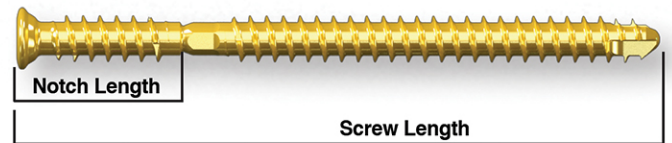

- **Lateral removal feature**
 - Hex feature allows for lateral removal of tibial fragment through the fibula

- **Medial removal feature**
 - Medial removal feature on screw tip allows for medial removal using a specialized driver

1. Internal data on file, TR-17102601

GBSS-01 RevA
™ Trademarks and ® Registered Trademarks of Paragon 28, Inc.
© Copyright 2017 Paragon 28, Inc. All rights reserved.
Patents: www.paragon28.com/index.php/patents

Paragon 28, Inc. 
4B Inverness Ct. E., Suite 280
Englewood, CO 80112 USA
(855) 786-2828

Paragon 28 Medical Devices Trading Limited
43 Fitzwilliam Square West
Dublin 2
D02 K792
Ireland
+353 (0) 1541 4756

For additional information on
Paragon 28® and its products
please visit: <http://www.paragon28.com>



- **Washer**
 - Available when the Breakaway screw is used without or outside of a plate

CE 0086

Exclusively foot & ankle
Paragon²⁸
®